



Now Available – NicVue v. 3.0.6

Patient Management Database

What is NicVue?

- A SQL-based multi-modality patient information and data management solution
- One database allows centralized management of all clinical neurophysiology modalities including Nicolet EEG, LTM, ICU, Sleep, EMG and IOM
- Supports legacy Nicolet data and databases allowing easy migrations
- Works with NicVue Connect to provide a bi-directional interface to any HL7-compliant hospital information system
- Role-based user management folders with audit trail for technicians, physicians and administrators
- Extensive patient demographic fields including 5 user-definable fields

v.3.0.6 Update

- Introduces NicVue Connect, a more reliable and less expensive HL7 solution
- Updates earlier versions of NicVue on devices such as Viking™, Synergy, NicoletOne™ and Endeavor™ CR
- Provides 64-bit operating system compatibility to support high technology workstations
- Easily convert reports to PDF locally and send to a repository for pick-up by your EMR
- Enhanced support for multiple orders per patient


Benefits

- In conjunction with NicVue Connect:
 - » Reduce transcription errors by importing demographics from your Hospital Information System
 - » Save valuable time importing rather than manually entering patient demographics
 - » Centralize clinical reports in your Electronic Medical Records
- Improve continuity of care by being able to see all clinical neurophysiology modalities in one place
- Improve staff productivity with the ability to manage all exam data from any NicVue workstation
- Save time on data management with one-button archiving to removable or online media for long-term data retention
- Easily print your own media labels
- Enhance your workflow when finding patients or data of interest using our date and name filters or create your own customized filters
- Ability to import patient images when required by Standard of Care
- Remotely review data over the Internet using a Citrix server to make timely patient care-related decisions

The top 17 leading US hospitals as rated by US News and World Report¹ choose Nicolet Brand EMG products – Viking, Synergy or EDX.

NicVue Connect HL7 Capabilities for Hospital Information System Integration

- Automatically import accurate patient demographics from an Admission, Discharge, Transfer message (ADT)
- Receive Orders with patient demographics and scheduling information
- Receive demographics with HL7 Query
- Send reports in the native format of your Nicolet device or PDF copies of those reports from any NicVue workstation to your EMR
- Native formats are dependent on the modality and may include Microsoft® Word, HTML, and Text

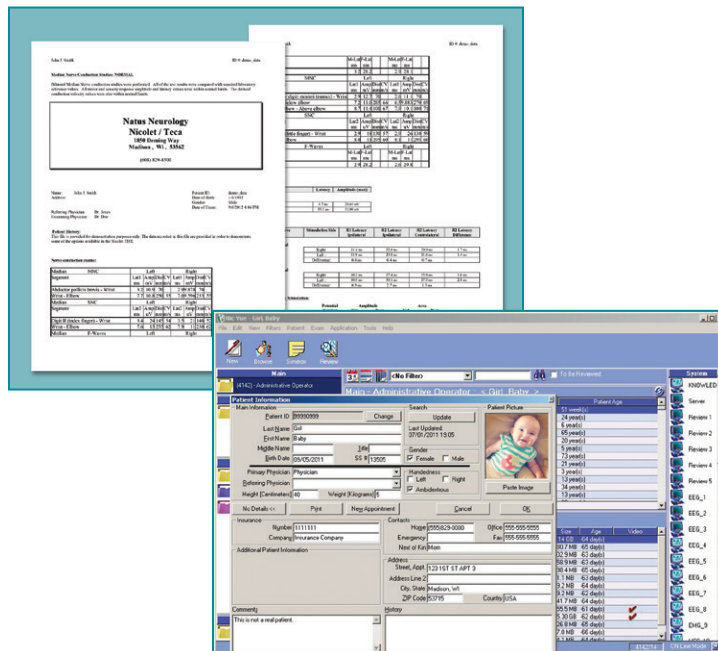

We have successfully implemented with the following systems:

- Epic®
- VistA®
- Cerner®
- Cloverleaf®
- McKesson®
- eGate™
- Sun®
- Siemens
- Patient.Net WCF

NicVue v.3.0.6 can be installed on new systems, updated from earlier versions of v.3.0, and upgraded from v.2.9 and earlier.

Nicolet Certified installation also available.

NicVue v.3.0.6 Recommended Software

System Platforms Required Application Software

NicoletOne	v.5.7 or newer
Viking	v.12 or newer
Synergy	v.15 or newer
EEG Vue	v.3.6 or newer
Endeavor CR	v.5.1.3 or newer

NicVue v.3.0.6 System Requirements:

Minimum

- Intel P4-3.2GHz
- 512MB RAM
- 750MB of free disk space on C:
- CD-R
- Windows® XP Professional (SP3)

Recommended

- C2-Duo 3+GHz
- 1024 MB RAM or better
- 1GB of free disk space on C:
- DVD-R
- Windows® 7

Server Requirements:

- Windows® Server 2003 or 2008 R2

Supplies

Natus Neurology offers a full range of neurodiagnostic accessories and supplies promoting patient comfort. Our dedicated customer service team provides a streamlined order and shipping process to save you time and money.

To learn more about Natus Neurology Service Programs or our full line of Supplies and Accessories, contact your local distributor or sales representative.
US Customers Call: 1-800-356-0007



Service

Natus Neurology is committed to providing exemplary service to our customers. Our dedicated and experienced Customer Service Team will assist with every aspect of an order. To support our products, we provide factory-trained Field Technicians and Clinical Application Specialists for onsite support. Additionally, we provide an in-house Technical Support Team, staffed with experts, and a strong distribution network in International Markets to offer a wide range of service options. Allowing our customers more time to care for their patients is our goal. Customer loyalty is our reward.



Natus Neurology Incorporated

3150 Pleasant View Road
 Middleton, WI 53562 USA
 Tel: 1-800-356-0007
 1-608-829-8500
 Fax: 1-608-829-8709
www.natus.com

1. Kimberly Leonard, "Best Hospitals 2013-14: Overview and Honor Roll", July 16, 2013, <http://health.usnews.com/health-news/best-hospitals/articles/2013/07/16/best-hospitals-2013-14-overview-and-honor-roll>

© 2014 Natus Medical Incorporated. All Rights Reserved. All product names appearing on this document are trademarks or registered trademarks owned, licensed to, promoted or distributed by Natus Medical Incorporated, its subsidiaries or affiliates.